

**A92117Mu02**

**Polyclonal Antibody to Activin A Receptor Type II A (ACVR2A)**

***Instruction manual***

FOR IN VITRO USE AND RESEARCH USE ONLY  
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

5th Edition (Revised in January, 2013)

**[ PRODUCT INFORMATION ]**

**Immunogen:** ACVR2A

**Clonality:** Polyclonal

**Host:** Rabbit

**Species Reactivity:** Mouse


**Ig type:** Rabbit IgG

**Purification:** Antibodies are purified by target protein affinity chromatography.

**Applications:** WB, ICC, IHC-P, IHC-F

**Form:** Liquid

**Size:** 100µg

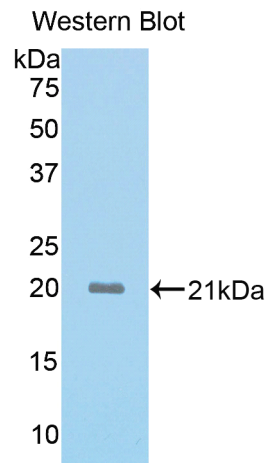

*Sample: Recombinant mouse ACVR2A*

**[ IMMUNOGEN INFORMATION ]**

**Immunogen:** Recombinant mouse ACVR2A (Lys326~Leu485) expressed in *E. coli*.

**Molecular Weight:** 19.6kDa

**USCN accession No.:** P92117Mu02

**Sequence:** The target protein is fused with N-terminal His-Tag and its sequence is listed below.

MGHHHHHSGSEF- KNVLL KNFLTACIAD FGLALKFEAG KSAGDTHGQV GTRRYMAPEV  
LEGAINFQRD AFLRIDMYAM GLVLWELASR CTAADGPVDE YMLPFEEEEIG QHPSLEDMQE  
VVVHKKKRPV LRDYWQKHAG MAMLCETIEE CWDHDAEARL SAGCVGERIT QMQRLL

### **[ ANTIBODY SPECIFICITY ]**

Anti ACVR2A is a rabbit polyclonal antibody raised against mouse ACVR2A. This antibody has been selected for its ability to recognize mouse ACVR2A in immunohistochemical staining and western blotting, non cross-reactive with other members of the family.

### **[ APPLICATIONS ]**

Western blotting: 1:100-400

Immunocytochemistry in formalin fixed cells: 1:100-500

Immunohistochemistry in formalin fixed frozen section: 1:100-500

Immunohistochemistry in paraffin section: 1:50-200

Optimal working dilutions must be determined by end user.

### **[ CONTENTS ]**

**Form & Buffer:** Supplied as solution form in PBS, pH 7.4, containing 0.02%Na<sub>3</sub>N, 50% glycerol.

### **[ QUALITY CONTROL ]**

**Content:** The quality control contains recombinant mouse ACVR2A (Lys326~Leu485) disposed in loading buffer.

**Usage:** 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemiluminescent (ECL).

Note: The quality control is specifically manufactured as the positive control. Not used for other purposes.

**Loading buffer:** 100mM Tris(PH 8.8), 2%SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, Na<sub>3</sub>N 0.02%.

### **[ STORAGE ]**

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.



Gentaur Molecular Products  
Voortstraat 49  
1910 Kampenhout, BELGIUM  
Tel 0032 16 58 90 45 | Fax 0032 16 50 90 45  
[www.gentaur-worldwide.com](http://www.gentaur-worldwide.com)  
[info@gentaur.com](mailto:info@gentaur.com)