



PAA605Hu01

Polyclonal Antibody to Adiponectin (ADP)

Organism: Homo sapiens (Human)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

8th Edition (Revised in Jun, 2013)

[PRODUCT INFORMATION]

Immunogen: ADP, Human

Clonality: Polyclonal

Host: Rabbit

Immunoglobulin Type: IgG

Purification: Affinity Chromatography.

Applications: WB, ICC, IHC-P, IHC-F, ELISA

Concentration: 200µg/mL

UOM: 100µg

[IMMUNOGEN INFORMATION]

Immunogen: Recombinant ADP (Tyr111~Asn244) expressed in *E.coli*.

USCN Accession No.: RPA605Hu01

Sequence: The target protein is fused with two N-terminal Tags, His-tag and S-tag and its sequence is listed below.

MHHHHHSSG LVPRGSGMKE TAAAKFERQH MDSPDLGTDD DDKAMADIGS EF-
YRSAFSVGL E TYVTIPNMPI RFTKIFYNQQ NHYDGSTGKF HCNIPGLYFF AYHITVYMKD
VKVSLFKKDK AMLFTYDQYQ ENNVDQASGS VLLHLEVG DQ VWLQVYGE GE
RNGLYADNDN DSTFTGFLLY HDTN

[ANTIBODY SPECIFICITY]

The antibody is a rabbit polyclonal antibody raised against ADP. It has been selected for its ability to recognize ADP in immunohistochemical staining and western blotting.

[APPLICATIONS]

Western blotting: 1:100-400


Immunocytochemistry in formalin fixed cells: 1:100-500

Immunohistochemistry in formalin fixed frozen section: 1:100-500


Immunohistochemistry in paraffin section: 1:50-200

Enzyme-linked Immunosorbent Assay: 1:100-200
Optimal working dilutions must be determined by end user.

[IMAGES]

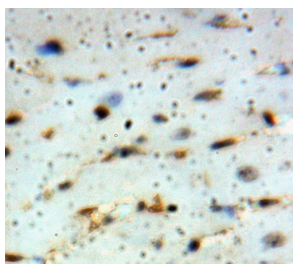

Used in DAB staining on formalin fixed
paraffin- embedded Heart tissue

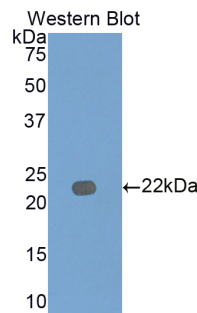

Used in Western Blot ,Sample:
Recombinant ADP,Human

[CONTENTS]

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN_3 , 50% glycerol.

[QUALITY CONTROL]

Content: The quality control contains recombinant ADP (Tyr111~Asn244) disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.
5uL per well when used in enhanced chemiluminescent (ECL).

Note: The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading Buffer: 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN_3 0.02%.

[STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.



GENTAUR

Gentaur Molecular Products

Voortstraat 49

1910 Kampenhout, BELGIUM

Tel 0032 16 58 90 45 | Fax 0032 16 50 90 45

www.gentaur-worldwide.com

info@gentaur.com