


**A91481Mu02**  
**Polyclonal Antibody to Cadherin, Neuronal (CDH2)**  
*Instruction manual*

FOR IN VITRO USE AND RESEARCH USE ONLY  
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

5th Edition (Revised in January, 2013)

**[ PRODUCT INFORMATION ]**

**Immunogen:** CDH2  
**Clonality:** Polyclonal  
**Host:** Rabbit  
**Species Reactivity:** Mouse  
**Ig type:** Rabbit IgG  
**Purification:** Antibodies are purified by target protein affinity chromatography.  
**Applications:** WB, ICC, IHC-P, IHC-F  
**Form:** Liquid  
**Size:** 100µg

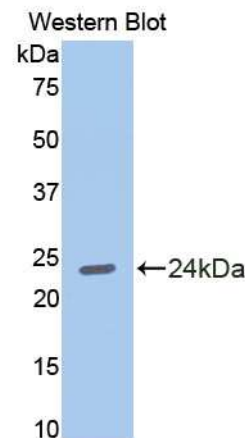

*Sample: Recombinant mouse CDH2*

**[ IMMUNOGEN INFORMATION ]**

**Immunogen:** Recombinant mouse CDH2 (Ala384~Glu562) expressed in *E.coli*.  
**Molecular Weight:** 21.3kDa  
**USCN accession No.:** P91481Mu02  
**Sequence:** The target protein is fused with N-terminal His-Tag and its sequence is listed below.  
MGHHHHHSGSEF-AMTFYGE    VPENRVDVIV    ANLTVTDKQ    PHTPAWNAAY  
RISGGDPTGR    FAILTDPNSN    DGLVTVVKPI    DFETNRMFVL    TVAAENQVPL    AKGIQHPPQS  
TATVSVTVID    VNENPYFAPN    PKIIRQEEGL    HAGTMLTTLT    AQDPDRYMQQ    NIRYTKLSDP  
ANWLKIDPVN    GQITIAVLVD    RE

Gentaur Molecular Products  
Voortstraat 49  
1910 Kampenhout, BELGIUM

Tel 0032 16 58 90 45 | Fax 0032 16 50 90 45  
www.gentaur-worldwide.com  
info@gentaur.com

**[ ANTIBODY SPECIFICITY ]**

Anti CDH2 is a rabbit polyclonal antibody raised against mouse CDH2. This antibody has been selected for its ability to recognize mouse CDH2 in immunohistochemical staining and western blotting, non cross-reactive with other members of the family.

**[ APPLICATIONS ]**

Western blotting: 1:100-400

Immunocytochemistry in formalin fixed cells: 1:100-500

Immunohistochemistry in formalin fixed frozen section: 1:100-500

Immunohistochemistry in paraffin section: 1:50-200

Optimal working dilutions must be determined by end user.

**[ CONTENTS ]**

**Form & Buffer:** Supplied as solution form in PBS, pH 7.4, containing 0.02%Na<sub>3</sub>N, 50% glycerol.

**[ QUALITY CONTROL ]**

**Content:** The quality control contains recombinant mouse CDH2 (Ala384~Glu562) disposed in loading buffer.

**Usage:** 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemiluminescent (ECL).

Note: The quality control is specifically manufactured as the positive control. Not used for other purposes.

**Loading buffer:** 100mM Tris(PH 8.8), 2%SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, Na<sub>3</sub>N 0.02%.

**[ STORAGE ]**

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.