



P91040Hu01

Vimentin (VIM)

Organism: Homo sapiens (Human)


Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY

NOT FOR USE IN DIAGNOSTIC OR THERAPEUTIC PROCEDURES

4th Edition (Revised in February, 2012)

[DESCRIPTION]

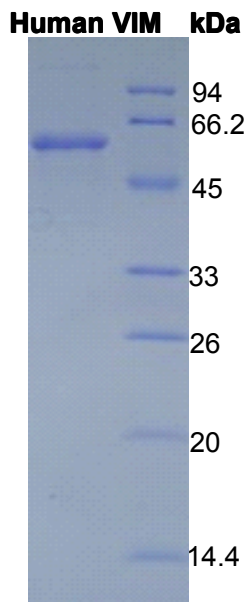

15% SDS-PAGE

Protein Names: Vimentin

Synonyms: VIM

Species: Human

Size: 10 μ g

Source: *Escherichia coli*-derived

Subcellular Location: Cytoplasm.

[PROPERTIES]

Residues: Ser2~Glu466 (Accession # P08670), with N-terminal His-Tag.

Grade & Purity: >92%, 55 kDa as determined by SDS-PAGE reducing conditions.

Formulation: Supplied as lyophilized form in PBS, pH 7.4, containing 0.01% Sarcosyl, 5% sucrose.

Endotoxin Level: <1.0 EU per 1 μ g (determined by the LAL method).

Applications: SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

Predicted Molecular Mass: 55.0 kDa

Predicted isoelectric point: 5.2

[PREPARATION]

Reconstitute in sterile PBS, pH7.2-pH7.4.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months .

Stability Test: The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed.(Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCES]

The target protein is fused with one N-terminal His-tag, its sequence is listed below.

MGHHHHHSGSEF-STRSVSSSS YRRMFGGPGT ASRPSSRSY VTTSTRTYSL GSALRPSTSR SLYASSPGGV YATRSSAVRL
RSSVPGVRL QDSVDFSLAD AINTEFKNTR TNEKVELQEL NDRFANYIDK VRFLEQQNKI LLAELEQLKG QGKSRLGDLY
EEEMRELRRQ VDQLTNDKAR VEVERDNLAE DIMRLREKLQ EEMLQREEAE NTLQSFRQDV DNASLARLDL ERKVESLQEE
IAFLKKLHEE EIQLQAQIQ EQHVQIDVDV SKPDLTAALR DVRQQYESVA AKNLQEAEEW YKSKFADLSE AANRNNDALR
QAKQESTEYR RQVQSLTCEV DALKGTNESL ERQMREMEEN FAVEAANYQD TIGRLQDEIQ NMKEEMARHL REYQDLLNVK
MALDIEIATY RKLLEGEESR ISLPLPNFSS LNLRETNLDS LPLVDTHSKR TLLIKTVETR DGQVINETSQ HHDDLE

[REFERENCES]

1. Brauksiepe B., *et al.* (2008) BMC Biochem. 9:25-25.
2. Honore B., *et al.* (1990) Nucleic Acids Res. 18:6692-6692.
3. Sommers C.L., *et al.* (1989) Cancer Res. 49:4258-4263.
4. Perreau J., *et al.* (1988) Gene 62:7-16.
5. Ferrari S., *et al.* (1986) Mol. Cell. Biol. 6:3614-3620.



Gentaur Molecular Products
Voortstraat 49
1910 Kampenhout, BELGIUM
Tel 0032 16 58 90 45 | Fax 0032 16 50 90 45
www.gentaur-worldwide.com
info@gentaur.com